

#### **About GE Healthcare**

GE Healthcare is a leading global medical technology and digital solutions innovator. GE Healthcare enables clinicians to make faster, more informed decisions through intelligent devices, data analytics, applications and services, supported by its Edison intelligence platform. With over 100 years of healthcare industry experience and around 50,000 employees globally, the company operates at the center of an ecosystem working toward precision health, digitizing healthcare, helping drive productivity and improve outcomes for patients, providers, health systems and researchers around the world.

Follow us on <u>Facebook, LinkedIn, Twitter and Insights</u>, or visit our website www.gehealthcare.com for more information.

- Work Related Musculoskeletal Disorders In Sonography, Society Of Diagnostic Medical Sonography, 2018, Susan Murphey, BS, RDMS, RDCS, CECD <a href="https://www.sdms.org/docs/default-source/Resources/work-related-musculoskeletal-disorders-in-sonography-white-paper.pdf?sfvrsn=8">https://www.sdms.org/docs/default-source/Resources/work-related-musculoskeletal-disorders-in-sonography-white-paper.pdf?sfvrsn=8</a>
- Forecasting the Future of Cardiovascular Disease in the United State, AHA Policy Statement, 2011, source: CIR.0b013e31820a55f5
- Source: Healthcare Infrastructure and Procedural Volume for Ultrasound Imaging, Frost & Sullivan, 2018. Approx. 108.12 million
  echo exams are performed annually. Calculation based on 26% total global prevalence of CVD cases (422 million) undergoing
  echo exam; extrapolated from US study indicating roughly 26% of total prevalent CVD cases underwent echo exams percentage
  value validated from reports.https://www.prb.org/wp-content/uploads/2015/12/2015-world-population-data-sheet\_eng.pdf
- 4. Kurt M, Shaikh K, Peterson L, et al. Impact on contrast echocardiography on evaluation of ventricular function & clinical management in a large prospective cohort. J Am Coll Cardiol. 2009; 53(9):802-810
- Depending on operation modes used and battery status
- 6. The Role of Al in Streamlining Echocardiography Quantification White Paper, Kristin McLeod JB80498XX
- Based on results of time and motion study conducted by GE "JB49055XX Cardiac Auto Doppler"; study results indicated time savings related productivity increase up to -8 on an annual basis for a facility per sonographer
- 8. European Association of Echocardiography recommendations for standardization of performance, digital storage and reporting of echocardiographic studies (Eur Journal of Echo 2008 Evangelista, Badano, Monaghan, Zamorano, Lancellotti).
- Recommendations for Quantification of Doppler Echocardiography: A Report From the Doppler Quantification Task Force of the Nomenclature and Standards Committee of the American Society of Echocardiography (JASE 2002)
- 10. Centricity Cardio Workflow v7 Intelligent Reporting out-of-the-box configuration compared to 2017 IAC guidelines excluding doppler, CCW Intellgent Reporting Outcome JB74831XX
- 11. EchoPAC Suite is a marketing name for EchoPAC Plug-in
- With the DICOM SR support, Measures & Analysis (M&A) for an exam can be sent at the end of the exam or when exported
  from local archive. The destination can be either a server on the network (Storage SCP) or a removable media (DICOM Media)
  depending on the DICOM dataflow selected. Custom measurements supported only for Adult Echo (TID5200) and Pediatric
  Heart (TID5220).
- 13. https://www.ncbi.nlm.nih.gov/pubmed/27689562
- $14. \quad \text{Data Breaches Will Cost Healthcare \$4B in 2019, Threats Outpace Tech, healthitsecurity, source: t.ly/xrAA}$
- 15. 5 Tips for Controlling Costs in Hospitals and Biomed Shops, source: t.ly/l9n7
- 16. Errors in Sonography, DOI: 10.1007/978-88-470-2339-0\_8.
- 17. e. a. M. Mårtensson, «High incidence of defective ultrasound transducers in use in routine clinical practice», European Journal of
- Echocardiography, vol. 10, no. 3, pp. 389-394, 2009. https://academic.oup.com/ehjcimaging/article/10/3/389/23966 https://probehunter.com/wp-content/uploads/FULLTEXT01.pdf
- 18. A multicentre survey of the condition of ultrasound probes, Ultrasound. 2016 Nov, Published online 2016 Aug 1. doi: 10.1177/1742271X16662301. https://www.ncbi.nlm.nih.gov/pmc/articles/PMC5098704/
- 19. GE internal data
- 20. Digital Expert is only offered in the USA

#### ©2020 General Electric Company. All rights reserved.


GE, the GE Monogram, Vivid, cSound, XDclear, HDlive, EchoPAC, ViewPoint, Centricity, Edison, iCenter, iLINQ and InSite are trademarks of General Electric Company. GE Healthcare, a division of General Electric Company.

DICOM is the registered trademark of the National Electrical Manufacturers Association for its standards publications relating to digital communications of medical information.

Ultra Edition is not a product name, it refers to the 2020 release of the Vivid portfolio

Third party trademarks are the property of their respective owners.

JB72911XX

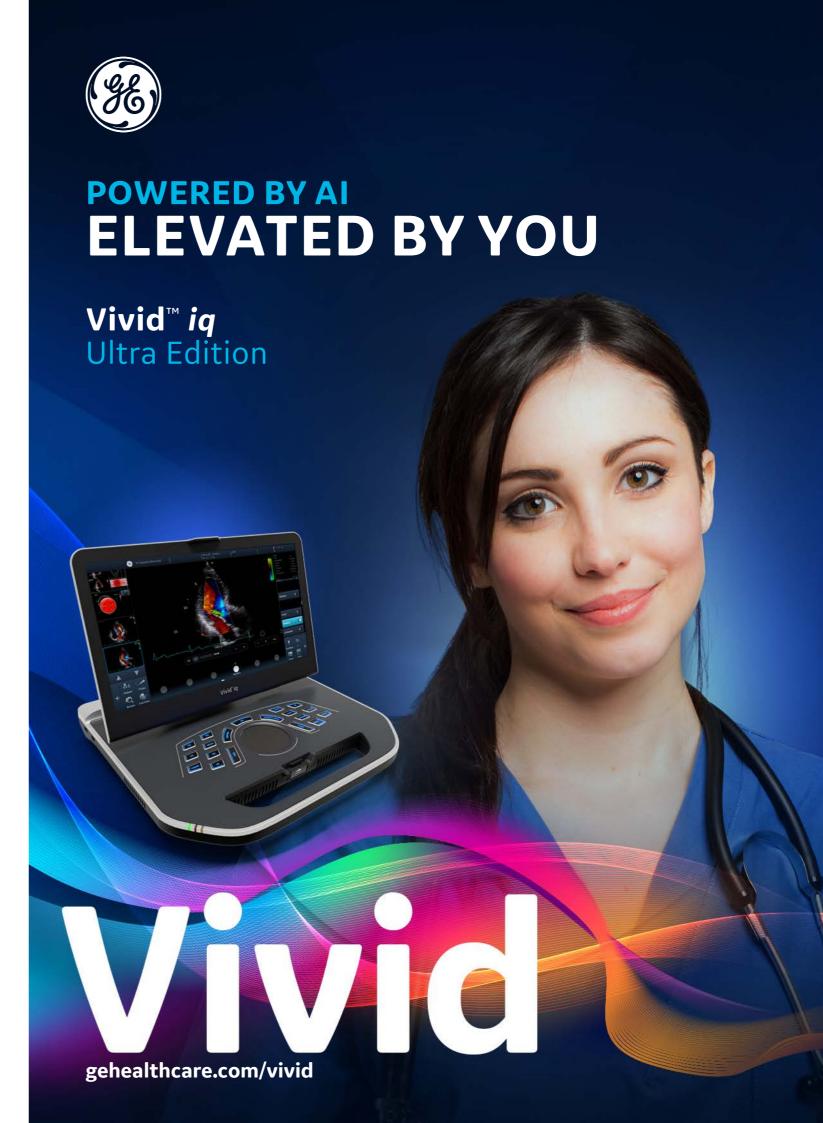


**490**%

of sonographers experience symptoms of Work Related Musculoskeletal Disorders (WRMSD)<sup>1</sup>

**\$120+** 

billion yearly

direct and indirect costs for employers<sup>1</sup>

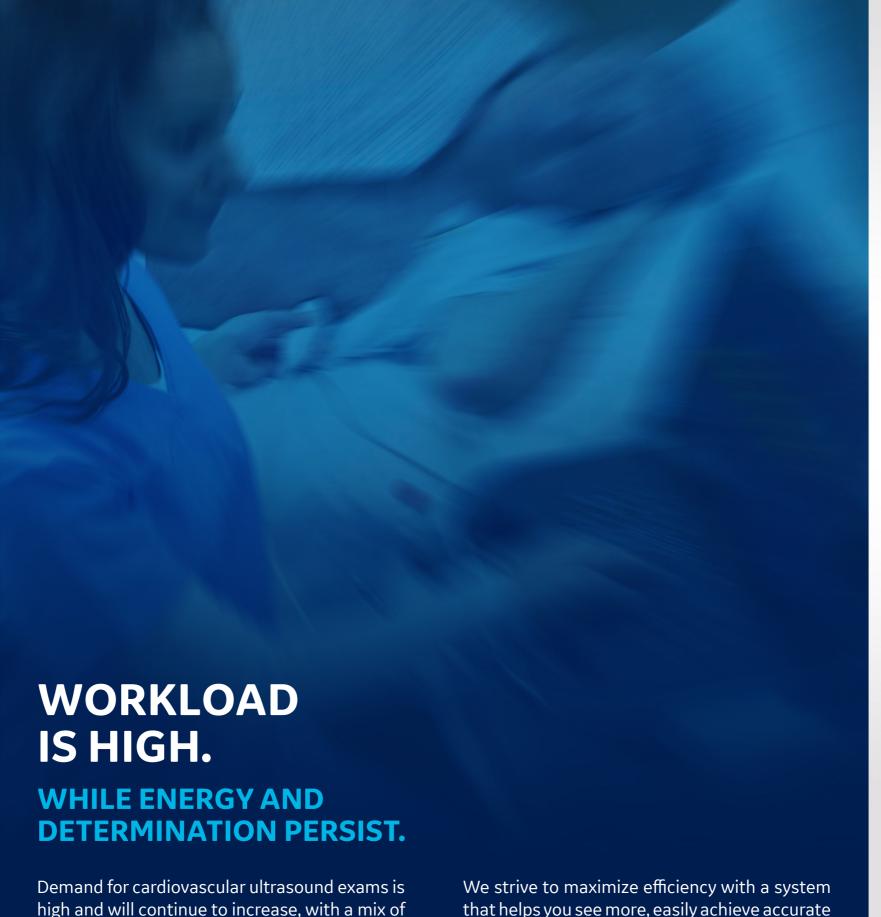

Using advanced clinical capabilities traditionally required extra effort and expertise - often resulting in delays and increased workload.

routine, follow-up and complex cases.

that helps you see more, easily achieve accurate measurements and minimize errors.

By 2030,

40.5%

of the US population is projected to have some form of Cardiovascular Disease (CVD)<sup>2</sup>

\_~108m

annual echo exams performed globally<sup>3</sup>

**10-15**%

of echo exams result in sub-optimal images<sup>4</sup>



### Best-in-Class Ergonomic Design

The improved ergonomic design is based on extensive user feedback and validation. The combination of the full touch screen control with new trackpad swipe and click gestures allow you to maintain support for your hand as you rest on a new ergonomic wrist rest. Scanning in comfort allows you to focus your attention on the image area.

Additional ergonomic features for your comfort include an adjustable monitor to customize viewing and typing angles, and a height-adjustable cart for comfortable scanning whether seated or standing.

### **Clinical Versatility**

Standard, Premium, 4D, and Point of Care versions are available, offering from basic 2D echo imaging to advanced 2D and 4D capabilities for quantifying heart and valve function, as well as ICE imaging and integrated connectivity for EP Lab. At the same time, diverse imaging presets and quantification tools for vascular and lung assessments are included, as well as a wide range of additional applications and presets.

### Ultra Portability and Flexibility

Designed to go places, the Vivid *iq* Ultra Edition is lightweight and has up to 4 hours of battery scan time on the cart. <sup>5</sup>

Even on the go, images can be processed anytime with the GE Healthcare Raw Data technology format. Once your scan is complete, share images quickly and safely between operators and readers with Wi-Fi connectivity and the Streaming Server feature.



90%

of sonographers experience symptoms of Work Related Musculoskeletal Disorders (WRMSD)<sup>1</sup>

**\$120**+

**billion yearly** direct and indirect costs for employers<sup>1</sup>





Do more on the trackpad



Minimize cognitive fatigue

# USER EXPERIENCE EXCELLENCE

### Your well being is our top priority

Prolonged scanning poses a known risk for work related musculoskeletal disorders (WRMSDs). Staff wellness is a key asset to your organization. Shorter exams, an intuitive user interface and adoption of an ergonomic body posture can positively impact your staff wellness. These are all possible with the Vivid *iq* Ultra Edition.

### A new level of user experience

The updated touch user interface was designed for convenience and fast scanning, for both left and right-hand users, minimizing muscle strain and helps sustain eye contact with the image area.

The new trackpad functionality and ergonomic wrist-rest accessory enable fast scanning and much needed ergonomic support.



Intuitive and familiar user interface

Your time is precious.
Save it.

DETECTABILITY<sup>6</sup> 98%

### AFI LV with AI View Recognition

Fully automatic recognition of the apical imaging views and measurements of GLS and segmental longitudinal Strain for LV.



## **POWERED BY AI**

### Improve diagnostic speed and accuracy

Vivid *iq* Ultra Edition features the latest AI-based technology to help reduce tedious tasks and improve workflow efficiency. You can diagnose more confidently and accelerate exams via automated (AI-driven) Cardiac Doppler and 2D LV measurements.

The results are impressive. Exam time is reduced, and operator fatigue minimized with up to 80% less clicks to get 2D measurements, and inter-observer variability diminished.<sup>6</sup>

Discover the many innovations brought by the Vivid *iq* Ultra Edition, and more importantly, how these can contribute to make clinical practice - Elevated by You.

Ultra Fast.
Ultra Precise.
Ultra Efficient.

### Al Cardiac Auto Doppler with Al

REDUCED TIME PER MEASUREMENT



Fewer Keystrokes<sup>7</sup>

LOWER INTER
OPERATOR VARIABILITY

**REDUCE VARIABILITY** 



Standardized exams with greater reproducibility <sup>7</sup>

ACCELERATED WORKFLOW



Productivity improvement

# **POWERED BY AI ELEVATED BY YOU** At GE Healthcare we strive to empower you by reducing wasted time and effort. We aim to remove tedious tasks and help make every moment count for your patients.

## CLINICAL EXCELLENCE for the Echo Lab

LESS CLICKS, UP TO 6




With the power of AI, the manual caliper measurements can be completed with 3 clicks: **Freeze - Measure - Auto.** A full set of reproducible measurements will instantly appear on the screen.



accuracy 6

### Al Auto Measure Spectrum Recognition

With the power of AI, a wide range of Doppler measurements can be completed with 2 clicks: **Freeze - Measure.** A Doppler trace and full set of associated measurements will instantly appear on the screen.

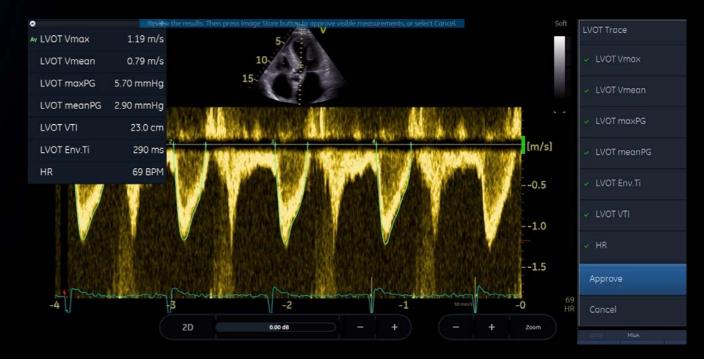



# SPEED AND VERSATILITY


### for Pediatrics


Speed and versatility are needed for your pediatric echo exams.

The Vivid *iq* Ultra Edition delivers high image quality with dedicated presets and advanced quantification tools, enabling precise measurements for all your patients – big or small.

### **Pediatric imaging**

Visualize small anatomies with speed, clarity and confidence thanks to Vivid *iq* Ultra Edition's superb high-resolution imaging and dedicated pediatric probes.

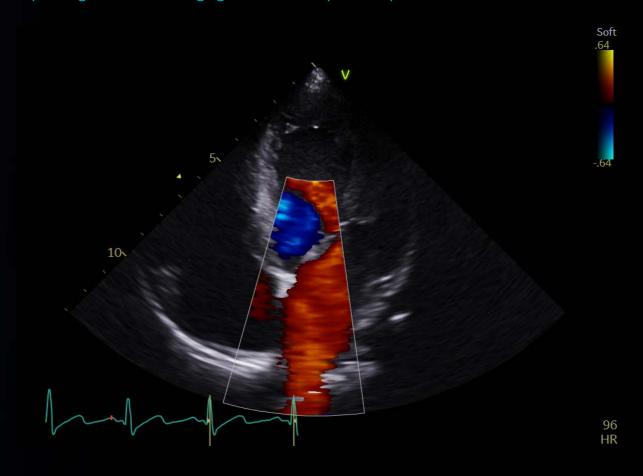

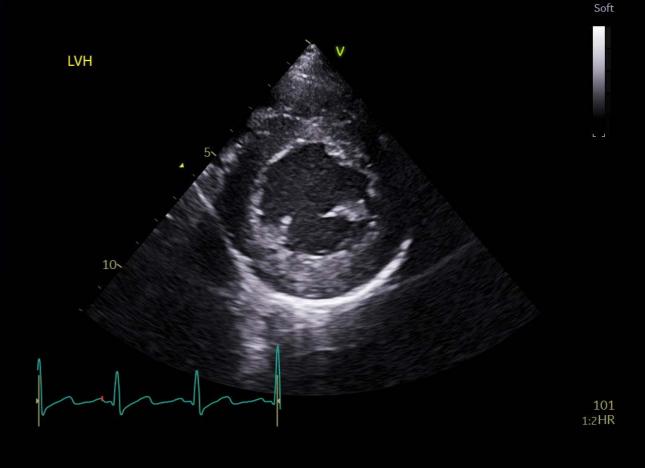



Use the Vivid *iq* Ultra Edition to assess heart and lung functionality in COVID-19 patients.


There are times when the best decision is to scan at the bedside to minimize risks. In addition to needing a portable system that is lightweight, battery powered, Wi-Fi enabled and easily cleanable, your demand for advanced imaging remains the same as in the Echo Lab. High-end performance for portable cardiac assessments is truly possible with the Vivid *iq* Ultra Edition.

With uncompromised heart visualization and excellent quantification tools, you'll save precious time with the GE Vivid AFI LV with Al-based View Recognition, and Auto EF for rapid LV Strain assessment. Al-powered measurements and tracing tools are easy to use, even by non-specialized users.

Use the new AFI RV tool to assess RV systolic function and estimate disease prognosis.


### **Lung Assessment**

Bedside lung assessment is made easier with dedicated FAST and lung presets that help to quickly optimize images.

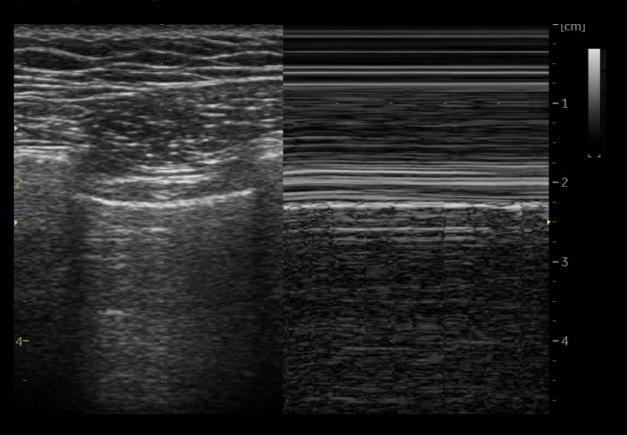

### **AFI RV**

Right Ventricle Global Strain (RVGLS) is an objective and accurate marker that is strongly correlated with RV ejection fraction. It has good prognostic power in many clinical studies.8

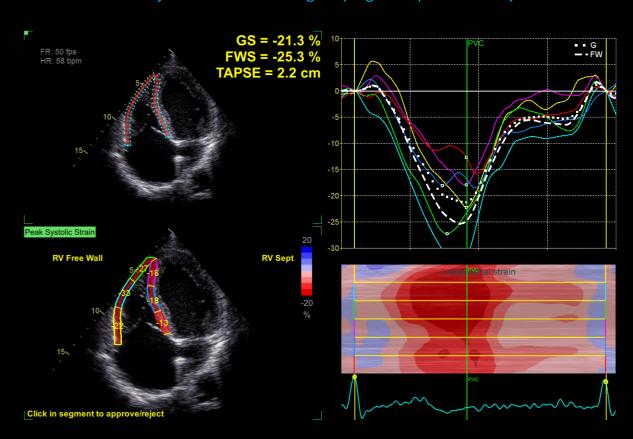



### **beyond Cardiology**


The demand for multi-purpose, cost efficient ultrasound systems with uncompromised image quality is growing. Your portable Vivid *iq* Ultra Edition will exceed your expectations across a wide range of applications.

With the aging population, vascular and abdominal exams are increasingly challenging. Vivid *iq* Ultra Edition provides you exceptional 2D and color flow imaging across applications including vascular, abdominal, fetal/obstetrics, and other applications beyond cardiology.

Efficient workflow, flexible touch screen display options (up/down, left/right), wide probe range and quantification tools ensure the confidence you need for speedy and accurate diagnosis also under challenging conditions.


### Vascular quantification

Develop fast and complete quantitative assessment of vascular anatomies, such as the Intima Media Thickness, with dedicated vascular measurement tools.

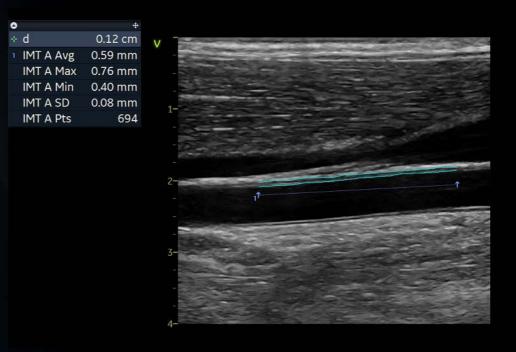

### **Abdomen diagnosis**

Visualize tissues and flow patterns with greater details thanks to Vivid iq Ultra Edition's high-resolution imaging.

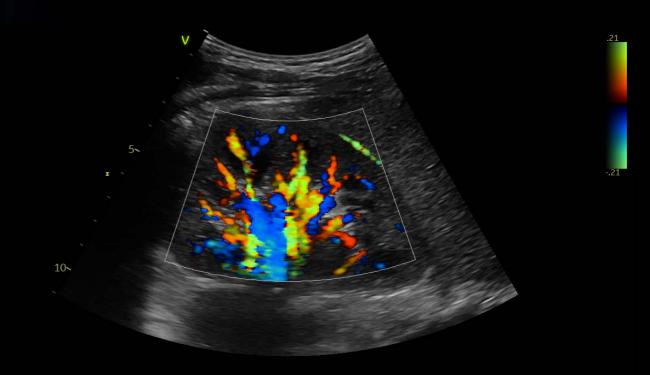



### for Interventional Procedures

When minimal footprint and infection control are critical, your heart team will appreciate a portable solution.

While intracardiac echo (ICE) is a dynamic and growing space, procedure rooms may be small and overcrowded, and infection control a constant struggle.


With the Vivid *iq* Ultra Edition minimal footprint, you can overcome your space constraints, while easily navigating imaging via the trackpad and touch screen, with your gloves on.

HD Color adds precision and insight to your everyday 4D TEE pre-operative screening and diagnostic procedures.

Expect more from your portable ultrasound and benefit from the Vivid family advanced multiplane and 4D visualization and navigation tools.

### **HD Color**

4D color flow rendering technique for semi-transparent visualization of origin and size of high velocity jets.

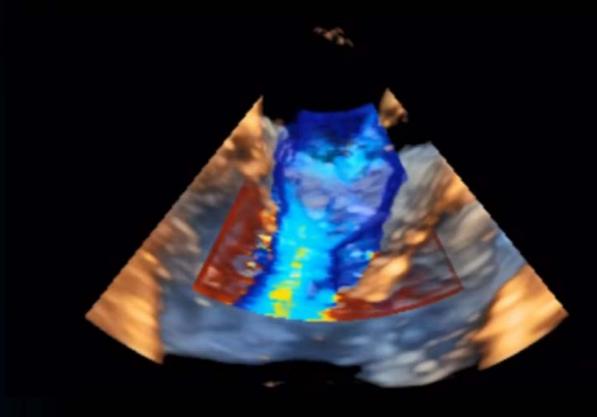

### 2D ICE

Great image quality, with easy and intuitive Intracardiac echo (ICE) touch user interface for your 2D ultrasound-guided ablation and structural heart procedures.



### Vivid iq Ultra Edition



### **VIVID HEART APPLICATIONS**

A wide range of clinical applications for use in Echo Lab, in the ICU and NICU, in the EP Lab and Cath Lab, and in Cardiovascular practices.




FLOW QUANTIFICATION Al Auto Measure Spectrum Recognition Al **VALVES AND CHAMBERS QUANTIFICATION** Al Auto 4D Auto AVO Measure 2D 4D Auto MVQ Auto EF 4D Auto LVQ

**FUNCTIONAL IMAGING** AFI LV with Al View Recognition Al **AFI RV** AFI LA



### VISUALIZATION AND NAVIGATION

### Why guess? When you can see.

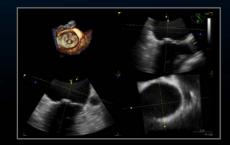


### **FlexiSlice**

With a distance gauge and two viewing layouts, this interactive tool for obtaining 2D or render views in live or replay mode may provide enhanced insight as well as save time.



### **FlexiViews**

Gain quick access to predefined 4D/ multiplane views during live mode, potentially reducing scan time during complex interventional procedures.

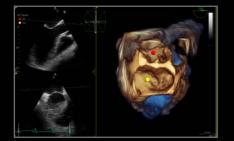

### **4D Markers**


Make annotations that are viewable from all angles on 4D ultrasound volume data sets and their 2D views, facilitating communication in the echo lab, cath lab and OR.



### View-X

See X-ray from fluoroscopy in real time right on your Vivid *iq* Ultra Edition screen as a picture in picture, facilitating communication between team members.

### New in Ultra Edition

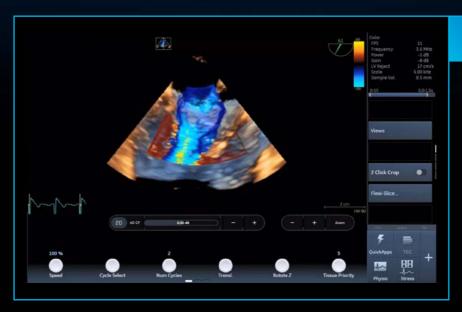


#### Ultra Edition

### **HD Color**

4D color flow rendering technique for semitransparent visualization of origin and size of high velocity jets.

#### Benefit:

- Enhance spatial relationships between flow and the surrounding structures
- Suppress non-diagnostic low flow information
- Work seamlessly with other visualization techniques such as 4D markers
- Supports 4D color flow data also from previous releases

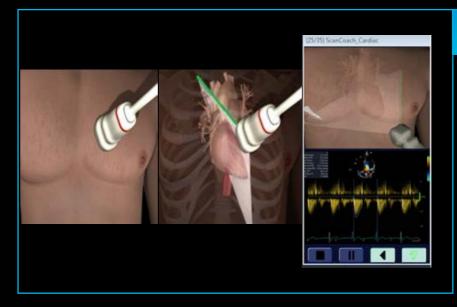

#### Ultra Edition

### **Scan Coach**

A reference and education tool that provides modules depicting basic scanning techniques with animated graphics of probe position, schematic of anatomy and reference clinical image.

#### Benefits:

#### **Guided scanning:**

- Conveniently located onboard the console within the Scan Assist Pro feature
- Provides guidance to obtain different views and measurements
- Helps inexperienced users performing echo scans


#### Education and reference tool:

- Assists in positioning the probe and probe orientation
- Anatomical reference helps visualize where the scan plane is located inside the heart
- Scan Assist Pro exam protocols can be customized per local guidelines and help ensure exam completeness



### **FLOW QUANTIFICATION**

Your time is precious. Save it.

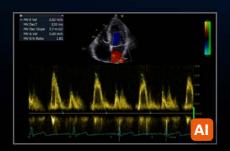


### **Cardiac Auto Doppler**

Semi-automatic Cardiac Doppler measurements.

#### Benefits

- Enhances reproducibility of follow-up studies when used in automated mode <sup>7</sup>
- Offers Doppler measurement in multiple cardiac cycles as recommended by guidelines for irregular heart rhythms <sup>8,9</sup>
- Supports less experienced users with advanced automation



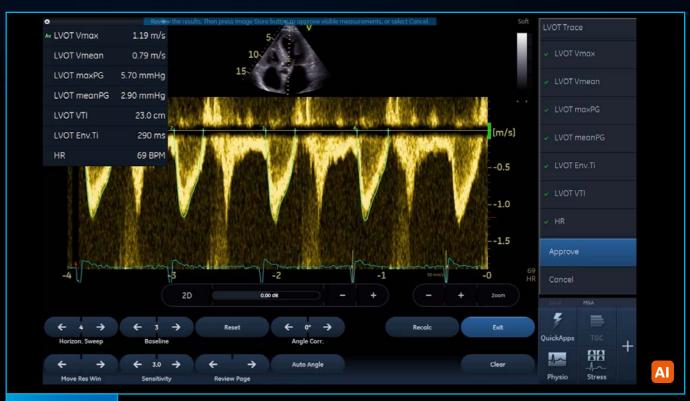

**Ultra Edition** 

### **Al Auto Measure - Spectrum Recognition**

Semi-automatic selection of appropriate spectral Doppler measurement tool.

#### Benefits:


- Enables fewer manual interactions by automatically opening the appropriate measurement tool  $^6$
- · Works seamlessly with Cardiac Auto Doppler
- Enhances reproducibility of follow-up studies when used in automated mode <sup>6</sup>
- Supports less experienced users with advanced automation



### **VALVES AND CHAMBERS QUANTIFICATION**


### New in Ultra Edition

### Precision at the heart of quantification.

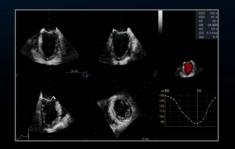


### **4D Auto LVQ**

Adapted to work with full volume data sets acquired with the 4D TEE transducer, 4D Auto LVQ for TEE brings you a fast, easy, two-click method of placing points to define the initial endocardial border.

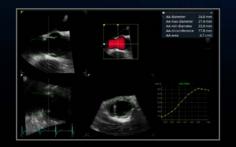


### **4D Auto AVQ**

Automatically segment, align and quantify the aortic outflow tract - vital to device sizing and orientation for TAVI/ TAVR procedures.

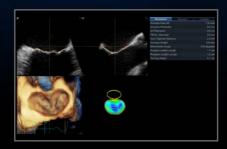

### **4D Auto MVQ**

Supporting TEE images, this integrated package helps visualize and quantify the mitral valve via a semi-automatic, surface-detecting algorithm.

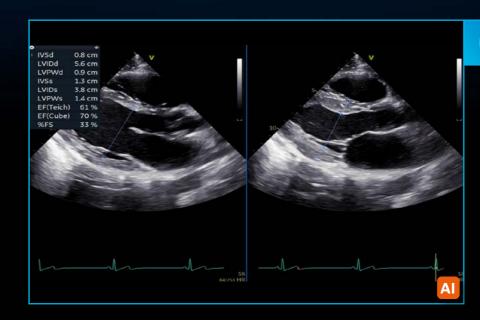

#### Ultra Edition

### **Al Auto Measure 2D**

Semi-automated LV dimension measurements (2D calipers) in the parasternal long axis view, reducing manual interactions.

- · Achieves fast measurements of left ventricle dimensions:
  - Up to 80% less clicks 6
  - No need to scroll to look for ED and ES frames
  - Reduce manual workflow during analysis of cardiac images
- Improves reliability and repeatability of measurements - potentially increasing reproducibility for follow-up studies





#### Ultra Edition

### **Auto EF**

Powered by Al-based View Recognition, Auto EF provides semi-automated quantification of left ventricular volumes and ejection fraction.

- Achieves fast measurements of ejection fraction
- DICOM® support. Assessment of the left ventricle ejection fraction also on data sets acquired on other vendors' systems




### **AFI FUNCTIONAL IMAGING**

From diagnosis to prognosis.




### New in Ultra Edition

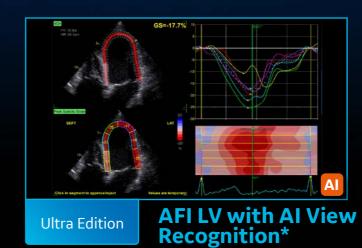

Powered by Al-based View Recognition, AFI LV provides semi-automated quantification of left ventricular global and segmental strain.

#### Benefits:


- Offers advanced industry pioneered speckle tracking algorithm for quantifying myocardial deformation
- Works seamlessly integrated ejection fraction calculation
- Supports Adult and Pediatric TTE and Adult TEE images
- Provides time savings via automatic selection of the appropriate 4-chamber, 2-chamber and APLAX images for analysis
- DICOM support. Assessment of the left ventricle ejection fraction also on data sets acquired on other vendors' systems

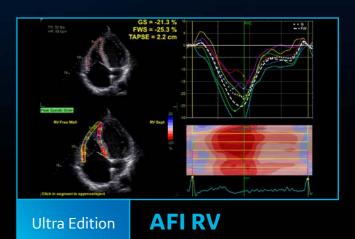

AFI RV is a novel tool to assess the right ventricular function by advanced speckle tracking echocardiography.

#### Benefits:

- Offers renown Vivid AFI user interface and workflow to allow current and new users easy adoption
- Supports right ventricle free wall strain, global strain and Tricuspid Annular Plane Systolic Excursion (TAPSE)
- Follows the 2018 EACVI-ASE Strain Standardized Task Force guidelines
- Supports right ventricle images also from previous releases

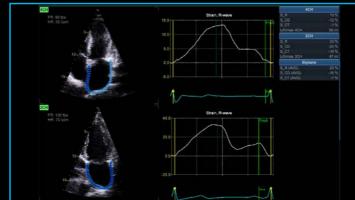


Ultra Edition

**AFILA** 

AFI LA Strain is a novel method to assess the left atrial function allowing global strain to be measured using speckle tracking echocardiography.

#### Benefits

- Offers Vivid renown AFI user interface and workflow allowing users to easily adopt
- Supports left atrium strain, volumes and emptying fraction measurements
- Follows the 2018 EACVI-ASE Strain Standardized Task Force guidelines
- Supports left atrium images also from previous releases

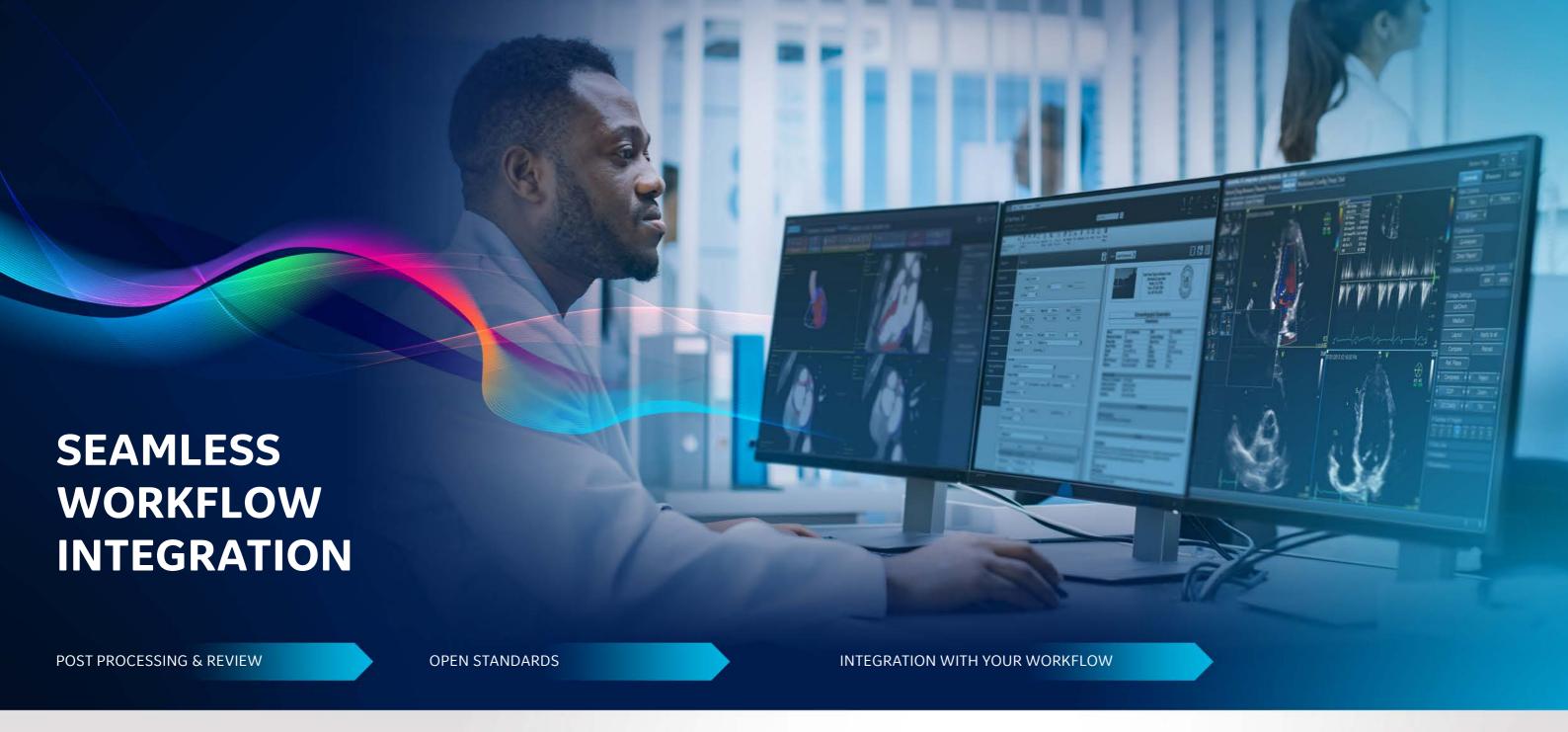

### **EchoPAC Software Only and EchoPAC Plug-in:**

- Analyze and review data from GE Healthcare Vivid family of scanners, as well as DICOM images from other ultrasound systems.
- Access all Vivid measurement and review tools utilizing GE Healthcare Raw Data or industry standard DICOM data
- DICOM Image transfer with optional GE Healthcare Raw Data transfers images easily in your existing workflow
- DICOM SR Measurement Transfer including standard and custom measurement allows seamless integration with GE Healthcare and other industry reporting systems and EMRs <sup>12</sup>

### EchoPAC Plug-in is available for:

- GE Healthcare Centricity™ Cardio Enterprise with Intelligent Reporting (IR).
- GE Healthcare ViewPoint<sup>™</sup> 6 with EchoPAC Suite <sup>11</sup>
- As a plug-in to third party PACS

With Centricity Cardio Enterprise IR, routine adult echo reports are

83%

complete before the physician opens the exam to review. 10



Protecting against these threats and safeguarding your patients and your institution requires more than anti-virus protection. SonoDefense is GE Healthcare's multi-layer strategic approach to cybersecurity and patient data privacy for ultrasound.

### SonoDefense is designed to:

- Keep the ultrasound machine safe and functional in the face of cyberthreats
- Protect patient data on the machine from unauthorized access
- Enable you to successfully implement patient data and security policies, while still managing product daily workflows

### SonoDefense strategy applies to Vivid portfolio:

- Windows® 10 IoT Secure Operating System provides multi-layered security
- Application whitelisting prevents malware execution
- Configurable user security provides user authentication and access control
- Data encryption protects stored data and during transmission
- Network firewall disables unneeded operating system services
- Integrates with existing facility security infrastructure



Healthcare is a soft-target for hacking and ransomware. 13

\$4B cost added in 2019 14



Ultrasound is especially vulnerable to operator-dependence leading to

### **VARIABILITY**

between exams 16



Constrained budgets increase pressure to

DO MORE with less 15

and to optimize assets



Probe mishandling can lead to damage which can cause faulty data and may lead to incorrect medical decisions 17,18



### [POP]

## PERFORMANCE OPTIMIZATION PARTNERSHIPS

Purchasing a GE Healthcare ultrasound is not only getting access to a high-technology or remarkable clinical applications. It's about enjoying a new user experience, at every step of ownership. We help you to outperform today, while preparing your department for tomorrow's challenges.

### STAFF EXCELLENCE

A comprehensive portfolio of training for clinical and technical users.

Helping you and your team build customized development plans to foster excellence and increased confidence.

### PROACTIVE MANAGEMENT

Use digital technology and tools to minimize expensive and disruptive unplanned downtime.

Proactive monitoring to help reduce cost and revenue loss from unplanned failures and automated updates for peace of mind.

### **DEVICE PROTECTION**

Keep your device state-of-the-art with software upgrades, new applications and security patches.

Optimizing your device to drive clinical and operational benefits and help you stay ahead of the game, without changing your equipment.

### **ASSET OPTIMIZATION**

Customizable dashboards for asset utilization and consulting services to provide actionable insights.

Achieving more with your assets to improve patient care and realize department strategic plans.

### **IMPROVED UPTIME**

Best-in-class repair services to drive uptime. Fully scalable from full coverage to shared maintenance.

Thoroughly aligned with your own in-house capabilities, providing the right balance between staff autonomy and our expertise.

### PROBE PERFORMANCE

Customizable portfolio of solutions for probe lifecycle needs to improve availability and performance.

Proactive probe care that may help you increase diagnostic quality, decrease cross-contamination risk and expand the life span of the transducers.

### Ready to make your Vivid POP?

Complete lifecycle solution for clinical, operational & financial outcomes. You take care of your patients, we'll take care of you.

### PARTNERSHIPS DRIVE RESULTS GE Healthcare is by your side to overcome these risks, helping you to:

- ✓ Keep your systems up and running, safe from breaches and cyberattacks
- Achieve more with your existing systems, without changing your investment plan
- ✓ Improve your activity, exam flow and staff planning, based on comprehensive data and reports
- Create comfort zone for your teams, reaching operational efficiency and clinical excellence
- Achieve high standard in probe-related cross-contamination and diagnostic errors

### **UNLEASH THE POWER** OF CONNECTED DEVICES

Your Vivid system has been designed to provide you with an optimal user experience. Connectivity is the key element to enjoy it fully, whenever and wherever you need it, regardless of site access restriction and planning constraints.

Discover a new world of services, included in every package:

### REMOTE TECHNICAL SUPPORT

Access to experts anytime, anywhere

**InSite™** connectivity enables OnDemand and realtime access to remote GE Healthcare experts

- Reduce disruptions
- Decrease system downtime
- Improve asset usage and staff productivity

It provides secure remote connectivity without requiring any open inbound ports or VPN connection.

#### **REMOTE FIX**

issues fixed remotely with InSite 19

### **PREDICTIVE MAINTENANCE**

Know the failure before it occurs

Transform unplanned downtime into planned service events with **OnWatch** technology. It provides automated, 24/7 system monitoring, capable of detecting a system failure before it occurs. Any deviation alerts our GE Healthcare engineers, who proactively work to keep your operations running smoothly.

#### **COST SAVINGS**





cost reduction, based on OnWatch alerts19

### **NEW REMOTE SUPPORT**

Real-time and interactive applications support

Digital Expert<sup>20</sup> and STAR provide an interactive, real-time, flexible & convenient way to get education and support.

- May help improve training outcomes
- Increase capacity and efficiency
- · Train staff on a short timeline



### **NEW AUTOMATED UPDATES**

No need to worry about your system safety

Get automated software updates with **eDelivery** along with safety patches enabled by remote software download. **SonoDefense** provides the highest level of

cybersecurity to keep your systems up-to-date, with no impact on your operations.

No on-site intervention needed.

### **DATA DRIVEN INSIGHTS**

All the insights you need to decide, at your fingertips

Better decisions start with better data **iCenter** is a secure, cloud-based asset management tool that offers comprehensive data analytics for your systems. It provides insights to make informed decisions and helps improve operational performance, optimize patient flow and maintain compliance standards.

**UpdateMe** is a complementary app that gives you access to the data 24/7 directly on your smartphone. You can receive notifications and create a service request anytime, anywhere.



